

THE ONE AND ONLY

Versatile, modular, and available in an endless number of configurations, the Quickie® 2 has been the industry's favorite ultra lightweight manual wheelchair for over three decades



QUICKIE 2

WITH NEW OPTIONS
FOR BOTH FUN
AND FUNCTION
THE QUICKIE 2
EMPOWERS YOU
WITH MORE
WAYS THAN
EVER TO
CONQUER
YOUR DAY



MODULAR FRAME DESIGN

The versatility of the Quickie 2 is in its modular frame that adapts with you as your function, body, and environment change. Independent, interchangeable frame components accommodate changing conditions without having to replace the entire wheelchair. You're always getting the best from your Quickie 2.




ENDLESS OPTIONS

With more options for rear wheels, casters, armrests, hangers, and footrests, and more opportunities to customize frame angles and measurements, the Quickie 2 is built to be an extension of you. Choose the hemi frame, depth-adjustable back, or angle-adjustable back for a more precise fit, or the Quickie Xtender™ Power Assist wheels for extra propulsion force.



NEW 5 Spoke Mag Wheels

This zero maintenance rear wheel option is a pound lighter per pair, features a closed spoke design to prevent debris build-up, and is available at no charge.

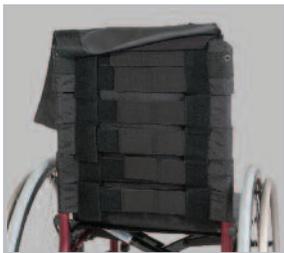


Tension Adjustable Back Upholstery

Choose from three tension adjustable upholstery options: 3DX, Ballistic, or the JAY® GO Back (pictured).



NEW Swing-In/ Swing-Out Hangers

The new hanger design features an intuitive bi-directional release lever and tapered locking mechanism to prevent loosening over time.

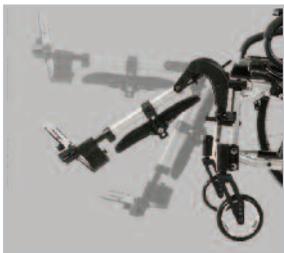

NEW Elevating/ Articulating Legrest

The improved height adjustable ELR/ALR allows for better alignment with the knee joint for a natural motion. Audible clicks index its position as it is elevated.





Q-FIT™ TECHNOLOGY

The Quickie 2 uses Q-Fit Technology to create a tight fit between components for a snug, out-of-the-box feel over the lifetime of the wheelchair. Living springs in the side guards and padded swing-away arms eliminate wiggle and rattling, while reversible seat saddles lock the seat rails in place for a smooth, rigid ride.



QUICKIE 2 LITE OPTION

Love the Quickie 2, but looking for something with an edge? The Quickie 2's minimalist counterpart, the Quickie 2 Lite, features an extremely lightweight, modular frame with a rounded rear for unique style and either fixed or swing-away front hangers.

* Conforms with ANSI/RESNA WC-19 standards. Please refer to Quickie® & Zippie® WC-19 Transit Securement Supplement (MK-100072) for additional information. Supplement is available for download at www.sunrisemedical.com

QUICKIE 2 SPECS

HCPCS Code: K0005

Overall Weight
29 lbs

Transport Weight
17.6 lbs

Weight Capacity
Standard: 265 lbs
Heavy Duty: 350 lbs

Seat Width
12" to 22"

Seat Depth
12" to 20"

Back Angle Adjustment
Fore (-) 3° to Aft (+) 15°

Front Seat Height
13.5" to 22"

Rear Seat Height
13" to 22"

Hanger Angle
60°, 70°, 80°, 90°

Transit Approved: Yes*

Warranty
Limited/Lifetime



©07.2014 Sunrise Medical (US) LLC
MK-100079 Rev.A

800.333.4000 www.SunriseMedical.com

