



Q6 edge[®] HD

With features like standard 4-pole motors, Mid-Wheel 6[®] Drive Design and ATX Suspension, the Q6 Edge[®] HD is engineered to meet the performance needs of the most active user.

The Q6 Edge HD accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value, and is available with our iLevel[®] technology.

quantumrehab.com

CA-Q6EDGEZ-87516



Beamsville, Ontario

888-570-1113

The Rehab Power Chair
Solutions Company

Specifications

DRIVE WHEELS	Standard: 14" Solid Optional: 14" Pneumatic	
CASTER WHEELS	Front: 6" Solid Rear: 6" Solid	
ANTI-TIP WHEELS	N/A	
MAXIMUM SPEED ⁴	Up to 4.7 mph	
GROUND CLEARANCE	3" (frame center)	
TURNING RADIUS ⁵	19.75"	
OVERALL LENGTH	35.6" without foot riggings	
BASE WIDTH	26.25"	
SEATING SIZES	TB2 HD	TB3 HD
Synergy[®]/Static Seat	Width: 20-28" Depth: 20-28"	Width: 20-32" Depth: 20-28"
Power Tilt	Width: 20-28" Depth: 20-28"	Width: 20-32" Depth: 20-25"
Lift & Tilt	N/A	Width: 20-24" Depth: 20-25"
SEAT-TO-FLOOR HEIGHTS	TB2 HD	TB3 HD
Synergy[®]/Static Seat ⁸	16.5-18.5"	17.5-19.5"
Power Tilt ⁸	18.2-20.9"	18.5-20.5"
Lift & Tilt ⁸	N/A	17.75-19.75"
MANUAL TILT/RECLINE ¹⁰	No	
BATTERY SIZE ⁶	Group 24	
BATTERY WEIGHT ¹⁴	52 lbs.	
AVAILABLE ELECTRONICS ¹²	75A Q-Logic 2 NE 75A Q-Logic 2 NE+ 100A Q-Logic 2 EX	
BATTERY CHARGER	8A Off-board	
MOTOR PACKAGES	4-Pole 4.7 mph ¹	
WEIGHT CAPACITY ⁷	450 lbs. ¹⁵	
BASE WEIGHT	143.5 lbs.	


- 1) Standard
- 2) Optional
- 3) Requires HD Suspension
- 4) Speed varies with user weight, terrain type, battery charge, battery condition and tire pressure.
- 5) Without front riggings
- 6) Battery capacity may change when adding power positioning systems.
- 7) Weight capacity is dictated by the selection of seating.
- 8) All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 300 lbs. weight capacity and solid tires (add 3/8" with pneumatic tires). Widths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.
- 9) 250 lbs. weight capacity with TRU-Balance[®] Systems
- 10) The Synergy[®] Manual Tilt weight capacity is 75 lbs. and maximum dimensions are 14" W x 14" D. Maximum dimensions of the Manual Recline vary by power base and seating system. Call for details.
- 11) Requires Group 34 batteries
- 12) All Quantum Rehab[®] units can be ordered with specialty controls
- 13) Required with high-speed motors
- 14) Battery weight may vary +/- 3 lbs. based on manufacturer.
- 15) Weight capacity is reduced to 400 lbs. with the iLevel[®] power adjustable seat height
- 16) A minimum seating size of 20"x20" is required with iLevel[®]. iLevel[®] maximum seating size is 24" x 24".

Features

- Mid-Wheel 6[®] allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance[®] 3 HD Power Positioning
- Low impact OMNI-Casters (nylon, spherical-shaped casters) on front and rear prevent wheel hang-ups
- ATX Suspension (Active-Trac[®] with extra stability) incorporates front OMNI-Casters and semi-independent rear caster beam for enhanced performance over more varied terrain
- iLevel[®] option enables up to 10" of seat elevation while at walking speed (up to 3.5 mph)^{15,16}

Options

- iLevel[®] power adjustable seat height^{15,16}
- TRU-Balance[®] 3 (450 lbs.) HD Tilt, HD Recline, and HD Tilt & Recline
- Q-Logic 2 Drive Control System
- USB mobile device charger
- Swing-away joystick mount
- New LED fender lights

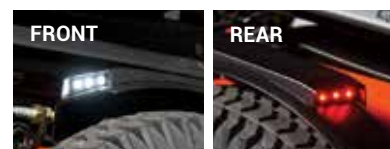

Colors

SEAT BACKS for units without TRU-Balance[®] 3 are only available in *Back in Black*.



Beamsville, Ontario • 888-570-1113 • www.quantumrehab.com • www.ilevel.rehab

