

www.hoggi.de

Made in Germany
BINGO
Evolution
by HOGGI 

BINGO was good - we just made it better



BINGO XPS
Extra Positioning System



BINGO IUS
Individual Upholstery System



BINGO Evolution - ASS
Adaptive Suspension System


BINGO Evolution

BINGO Evolution - Innovative Features

The **BINGO Evolution** is the result our goal of continuous product development. The **BINGO** was good - we just made it better!



Compact folding size

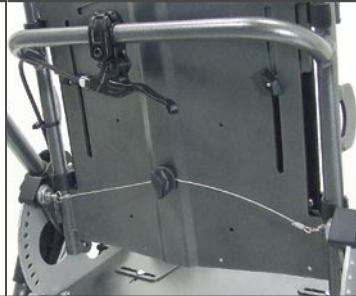

Tilt & backrest recline:
Central and easy to handle



Detachable seat



Seat in rear facing

BINGO Evolution - IUS & Frame & Upholstery colour



Frame colour: sparkling iron effect



Upholstery colour: Blue / Grey



Upholstery colour: Red / Grey

BINGO Evolution - XPS Seat - Extra Positioning System

XPS Features:

The XPS seat features the same great **Bingo Evolution** base with a seating system designed for moderate to involved seating needs.

- The shorter back allows headrests to be mounted in front or behind the plane of the backrest.
- Lateral pads provide aggressive positioning support and are adjustable from behind without transferring the child.
- Proper pelvic positioning achieved with the use of a contoured wheelchair cushion.



BINGO Evolution

BINGO Evolution - accessories



Bus Transport Package
Standard (ISO 7176-19)



Soft Accessory Bag
Standard (max. 3 kg)



Retainer standard with
grab rail or tray



Therapy tray grey or
clear IUS & XPS



Grab rail with upholstery
IUS & XPS



Five-point harness
IUS only



Groin strap
IUS & XPS



Chest-shoulder harness
IUS only



Abduction block
IUS only



Headrest pads
IUS only



Lap belt
IUS & XPS



Lateral trunk supports
IUS only

COBRA - Indoor base

Indoor base for **BINGO Evolution** IUS and XPS Seats. Features include:

- 45° of Tilt.
- Infinitely adjustable gas spring assist for height adjustment.
- Can be configured for the weight of the child.



BINGO Evolution

BINGO Evolution & COBRA - Technical data

Mobility base	Size 1		Size 2	
Push handle height	75-110 cm	29.5-43.3"	75-110 cm	29.5-43.3"
Overall width	63 cm	24.8"	68 cm	26.8"
Overall length mobility base	92 cm	36.2"	92 cm	36.2"
Wheel diameter front	19 cm	7.5"	19 cm	7.5"
Wheel diameter rear	29 cm	11.4"	29 cm	11.4"
Tilt in space	45° in front or rear facing		45° in front or rear facing	
folded size with rear wheels	80x63x38 cm	31.5x24.8x15"	80x68x38"	31.5x26.8x15"
folded size without rear wheels	80x63x31 cm	31.5x24.8x12.2"	80x68x31"	31.5x26.8x12.2"
Weight mobility base	8.6 kg	19 lbs.	8.7 kg	19.2 lbs.
Load capacity mobility base	60 kg	132 lbs.	60 kg	132 lbs
IUS Seat	Size 1		Size 2	
Seat Depth	16-34 cm	6.3-13.4"	22-40 cm	8.7-15.7"
Seat Width w/ hip support	18-30 cm	7.1-11.8"	23-35 cm	9.1-13.8"
Lower Leg Length	19-37 cm	7.5-14.6"	19-40 cm	7.5-15.7"
Back height	45-63 cm	17.7-24.8"	54-75 cm	21.3-29.5"
Recline	90° - 160°		90° - 160°	
Maximum load	50 kg	110 lbs.	50 kg	110 lbs.
Weight	6.4 kg	14 lbs.	7.4 kg	16.3 lbs.
XPS Seat	Size 1		Size 2	
Seat Depth	16-34 cm	6.3-13.4"	22-40 cm	8.7-15.7"
Seat Width w/ hip support	18-30 cm	7.1-11.8"	23-35 cm	9.1-13.8"
Lower Leg Length	19-37 cm	7.5-14.6"	19-40 cm	7.5-15.7"
Back height	36-52 cm	14.2-20.5"	41-60 cm	16.1-23.6"
Recline	90° - 160°		90° - 160°	
Maximum load	50 kg	110 lbs.	50 kg	110 lbs.
Weight	6.4 kg	14 lbs.	7.4 kg	16.3 lbs.
COBRA				
Seat Height	27 - 61 cm		10.6 - 24"	
Total width	61 cm		24"	
Total length	71 cm		28"	
Caster diameter	75 mm		3"	
Tilt	0° to 45°			
Total weight	11.3 kg.		25 lbs.	
Load capacity	50 kg.		110 lbs.	



Proudly distributed in Canada by:

R&M Health Care Inc. • 790 Redwood Square • Unit 8 , Oakville, ON L6L 6N3

Ph: 905-847-1983 • Toll free: 866-664-1983 • Fax: 905-847-2708

www.RandMHealthCare.ca