

permobil

THE POWER OF MOBILITY



M300
CORPUS HD


TOUGH ENOUGH



AVAILABLE FABRIC



LEATHERETTE

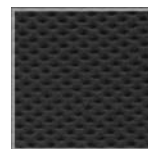

MESH

M300

CORPUS **HD**

AVAILABLE SHROUD COLORS



REALTREE AP CAMO



ONYX BLACK



VOLCANO RED



SUBLIME GREEN



PLATINUM SILVER



COBALT BLUE



PURPLE RAIN



PHANTOM GREY



POP STAR PINK

SPECIFICATIONS

Maximum User Weight	450 lbs	205 kg	Seat to Floor Height	(Elevator)	18", 19", 20"	460, 485, 510 mm
Maximum Speed	5 MPH	8 km/h		(Seat Tube)	18", 19", 20"	460, 485, 510 mm
Driving Range ¹	17.8 miles	28.6 km	Tilt Options	(Power)		0° - 45°
Base Width	26"	660 mm		(Manual)		No manual option
Base Length	38"	965 mm	Recline Options	(Power)		85° - 150°
Minimum Turning Radius	22"	555 mm		(Manual)		85° - 120°
Ground Clearance	3"	75 mm	Legrest Elevation Options	(Power)		85° - 175°
Optional Seat Elevator	8"	200mm		(Manual)		85° - 175°
Battery Type	Group 24, Group 34		Armrest Pad Lengths		13", 16", 18"	335, 410, 460 mm
Drive Electronics	120A R-net		Seat Widths		19", 21", 23", 25"	470, 520, 570, 620 mm
Maximum Slope ²	7.5° Up / 7.5° Down		Seat Depths		19" - 25"	470 - 620 mm
			Backrest Heights		21", 24" - 29"	530, 600 - 730 mm

¹ Actual driving range will vary based on driving and battery conditions.

² Refer to Owner Manual for safe driving practices



TOUGH ENOUGH



INTRODUCING THE NEW M300 CORPUS® HD, a heavy duty power wheelchair designed specifically for users up to 450 lbs. The HD has extremely sturdy arm and legrests as well as wider footplates for additional room, and it features the world-renowned Corpus seating system for the ultimate in comfort with an unmatched power tilt range of 0-45° and an unprecedented recline of 85-150°. The M300 base is highly maneuverable and has been both extended and widened for the HD model, while the suspension has been strengthened for an overall more stable platform. The ultra-strong actuators allow users to safely elevate their feet above their hearts, while the durable mesh-style upholstery helps keep moisture away from the body. With the legendary power and toughness of Permobil, the M300 Corpus HD is **BUILT TO WITHSTAND THE RIGORS OF ANY HEAVY DUTY LIFESTYLE.**

www.permobil.com




Power tilt, recline and elevating legrests allow users' feet to be positioned above the heart.



EVERYTHING YOU'VE COME TO EXPECT FROM THE CORPUS® SEATING SYSTEM. ONLY STRONGER.

0-45° Power Tilt
85-150° Power Recline
Power Elevating Legrests

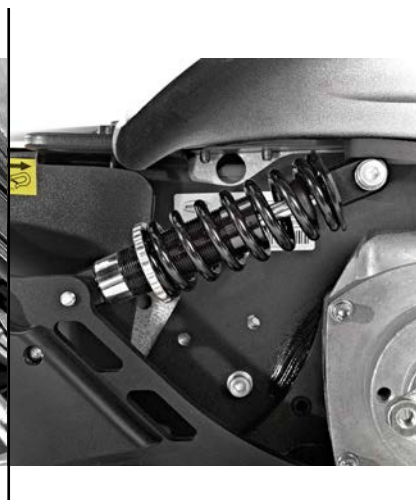

Adjustable legrests

Extremely durable, sturdy legrests

Adjustable footplates

Oversized, independently adjustable footplates

Strong actuators

Incredibly strong actuators for seat functions

Suspension

Independently adjustable suspension for a smooth ride under any condition