

QUICKIE


SR45TM

A well appointed economical option for tilt in space.

Proven durability

Proven functionality

Pre-configured for ease of ordering

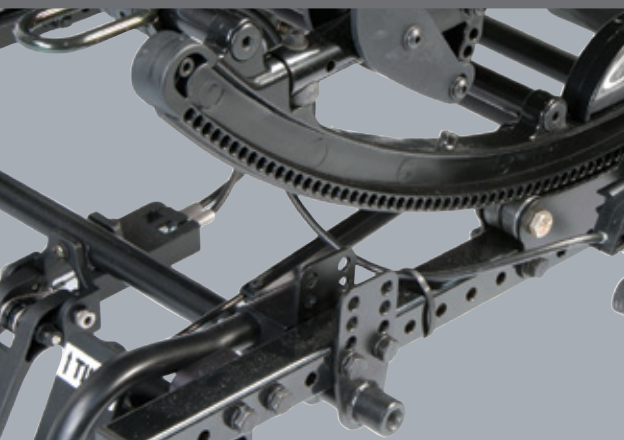

Nothing Beats a Quickie™

The Quickie® SR45™ is the newest member of the acclaimed Quickie tilt-in-space family and features 45 degrees of patented and proven IRIS tilt technology and a targeted offering of clinically effective options. Available in 3 pre-configured sizes straight off the shelf. The SR45 combines the best mix of options in a ready to go economical package.



Rotation-in-Space Technology

The Quickie SR45's seat frame rotates around the client's centre of gravity providing the shortest possible wheelbase without compromising stability. When the centre of gravity and centre of rotation are properly aligned, tilting is smooth, easy and controlled.

Smart Option Mix

SR45 comes standard with:

- Non folding backrest
- Height adjustable back with adjustable stroller handles
- Foot tilt
- Transit tie downs
- Dual post height adjustable flip back with full length arm pads
- 70` Swing in/ out hangers, solid seat pan and rear anti tips
- Composite foot plates
- 12" Mag rear wheels

These features will meet the needs of a wide range of clients. Available ex-stock as a quick ship model.

Newly designed for SR45



Low Profile Dual Post Armrest



Swing-in/Swing-out Hanger



Steel Base



Foot Release Tilt Actuator



SR45 with additional seating

SPECIFICATIONS Quickie SR45		QUICK SHIP SKU'S	
Frame style	Tilt-in-Space	SR45 16	16x16
Tilt range	0" – 45"	SR45 18	18x18
Product weight	18 kg	SR45 20	20x18
User weight capacity*	136 kg		
Frame colour	Black	Warranty	Lifetime on frame, 1 yr on parts
TGA #	100167	Transit approved	WC 19

* Please note: This includes user weight and items on chair.