




Features, Performance and Travel Ready

The Go-Go® Sport delivers high-performance operation and easy disassembly for convenience on-the-go. With a 325 lbs. weight capacity, charger port in tiller, and standard front LED lighting, the Go-Go Sport is feature rich and travel ready.







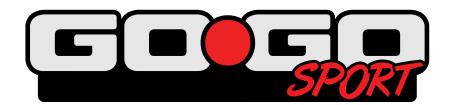

SPECIFICATIONS

Model number	S73 (3-Wheel)	
Model number	S74 (4-Wheel)	
Max. weight capacity	325 lbs.	
Maximum speed ²	Up to 4 mph	
Ground clearance ¹	3"	
Turning radius ¹	36" (3-Wheel)	52" (4-Wheel)
Overall length ¹	43" (3-Wheel)	43" (4-Wheel)
Overall width ¹	20"	
TIRES:		
Front	3W: 3" x 8" solid 4W: 3" x 9" solid	
Rear	3" x 9" solid	
A 11 11	_	

Anti-tips Rear

Standard seating Foldable, compact, vinyl-covered

molded plastic

Drivetrain 24-volt DC motor,

sealed transaxle

Dual braking system Regenerative and

electromechanical

Battery charger Off-board, 2 amp

Per-charge range^{2,5} 3-Wheel: Up to 9 miles

4-Wheel: Up to 12 miles

Battery requirements³ (2) 12V deep-cycle (standard)

Compatible batteries 18 AH

WEIGHT:

Batteries (2 included) 4 18 AH (27.4 lbs. battery pack)

Weight w/o battery pack 3-Wheel: 83.2 lbs.

4-Wheel: 99.2 lbs.

Heaviest piece 3-Wheel: 36.2 lbs. (front section)

4-Wheel: 52.2 lbs. (front section)

Interchangeable shroud panels





Red Blue


FEATURES


- Feather-touch disassembly permits simple frame separation with only one hand
- · Pride's exclusive black, non-scuffing tires
- Auto-connecting front-to-rear cable
- Includes 2 sets of easily changeable colored panels in Red and Blue
- Front frame-mounted seat post offers maximum stability
- All-new frame design easily disassembles into 5 super lightweight pieces for convenient transport and storage
- Modular design for easy serviceability
- Microprocessor-based controller offers optimal power management and added safety features
- Standard front tiller-mounted basket included
- Deck basket for additional storage (4-Wheel only)
- Convenient off-board dual voltage charger can charge battery pack on-board or off-board
- Charger port in tiller

Pride Lifts & Ramps

These Pride® Lifts are compatible with the Go-Go® Sport.

- Backpacker® Plus
- Backpacker® MV
- Power Tote Mini (3-wheel only)
- PrideBoom 250

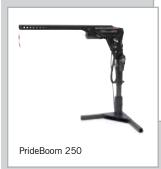

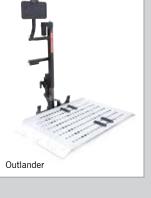

¹Due to manufacturing tolerances and continued product improvement, this specification is subject to a variance of +/- 3%.
²Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition.

This conditions are because the product of th

This specification can be subject to a variance of (+ or -) 10%. 3AGM or gel cell type recommended.

Battery weight may vary based on manufacturer.

STested in accordance with ANSI/RENSA, WC Vol 2, Section 4, and ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight canacity.

The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.



