

# PRISM TRUEFIT BACK




The Prism TrueFit™ was designed for users with a variety of positioning needs such as Kyphosis, Lordosis, Scoliosis or protruding gibbous deformities. The ability to customize the back to each individual's needs, plus the exclusive Prism hardware, with ease of adjustment and installation, contributes to the unique and personal "true fit" in our TrueFit Back!

Prism Series Backs		STYLE: TrueFit Back	
Stock Sizes – Same Day Shipment			
ADP CODE: SEBFM0050			
Part Number	Size W x H	Part Number	Size W x H
TFB1616	16" x 16"	TFB2016	20" x 16"
TFB1816	18" x 16"		
Non-Stock Sizes – Shipment within 3 business days			
ADP CODE: SEBFM0055			
Part Number	Size W x H	Part Number	Size W x H
TFB1618	16" x 18"	TFB1818	18" x 18"
TFB1620	16" x 20"	TFB1820	18" x 20"
TFB1716	17" x 16"	TFB2018	20" x 18"
TFB1720	17" x 20"	TFB2020	20" x 20"
Custom Sizes - Shipment within 7 business days			
ADP CODE: SEBFM0065			
Part Number	Size W x H		
TFBCUS	Custom		



Prism Series Backs		STYLE: TrueFit Back	
Bariatric Sizes – Shipment within 3 business days			
ADP CODE: SEBFM0060			
Part Number	Size W x H	Part Number	Size W x H
TFB2116	21" x 16"	TFB2316	23" x 16"
TFB2118	21" x 18"	TFB2318	23" x 18"
TFB2120	21" x 20"	TFB2320	23" x 20"
TFB2216	22" x 16"	TFB2416	24" x 16"
TFB2218	22" x 18"	TFB2418	24" x 18"
TFB2220	22" x 20"	TFB2420	24" x 20"
Custom Bariatric Sizes - Shipment within 7 business days			
Part Number	Size W x H		
TFBBCUS	Custom		

## Product Highlights:

- ✓ ½" of comfort foam layers on a ½" of support foam covering "ice cube" pattern of standard 2" x 2" foam blocks
- ✓ Pins adjustable to accommodate up to 2" difference in seat width (standard, 0.5" and 1" longer)
- ✓ Back angle adjustment forward 15 degrees, back 20 degrees, optimizing the users' positioning
- ✓ Extra foam covering prevents any exposure of shell
- ✓ Mix & Match foam cubes by using 2" polyurethane foam and ordering our optional 3" polyurethane or 2" viscos foam cubes
- ✓ Lightweight, sturdy aluminum shell
- ✓ **Product Weight:** 5.5 lbs (7.2 lbs bariatric), (18"W x 16"H)
- ✓ **User Weight:** 250 lbs (350 lbs bariatric)
- ✓ **Warranty:** 3 years limited on back, 90 days on cover, and 1 year on hardware

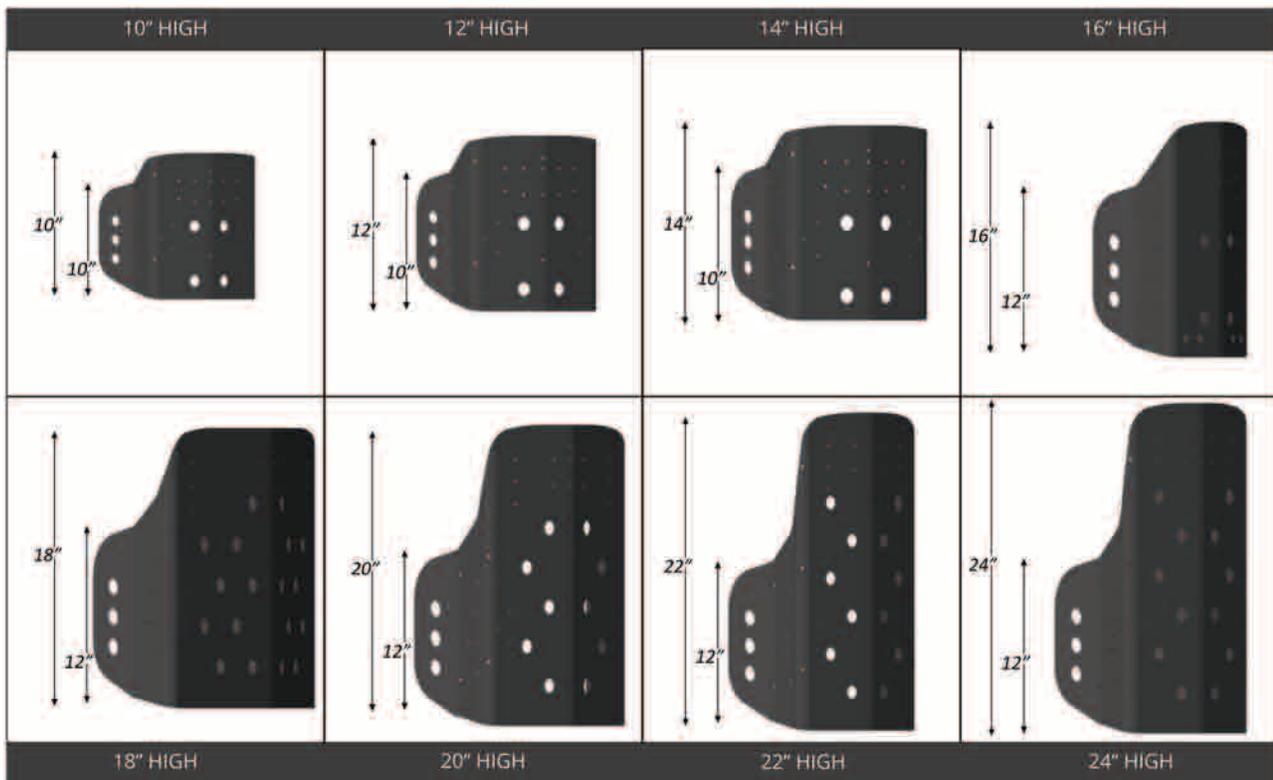

Lateral Dimensions - Truefitt Back		
TRUEFIT, Heights 12-24"		
Chair Width	DIM B (")	DIM A (")
14"	14	7.6
16"	15	7.6
18"	16	7.6
20"	17	7.6
22"	18	7.6

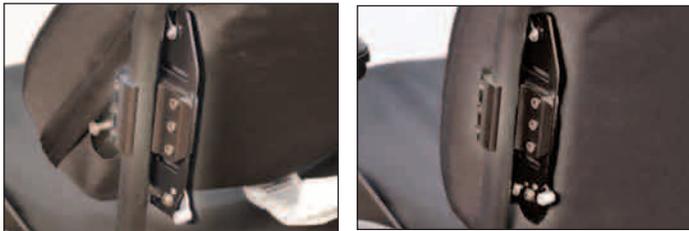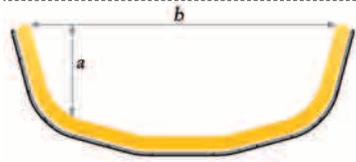

Cover, Foam Layers and Aluminum Hardware

- ✓ Fully adjustable for height, depth and angle of recline
- ✓ Hardware fits 3/4", 7/8" and 1" tubing and can be installed within a few minutes with the right tools
- ✓ Quick release lower latch levers allows the backrest to be easily removed.